

# EXCEPTIONAL PRODUCTIVITY IN A COMPACT DESIGN

A powerful benchtop system that maximizes valuable laboratory space

**Access 2 Immunoassay System**



» Move healthcare forward.

# » **PROVEN ANALYTICAL PERFORMANCE. SPACE-SAVING VALUE.**

**The Access 2** contains the robustness of a reference laboratory immunoassay analyzer in the convenient size of a benchtop system, providing **greater control over performance, reliability and speed without sacrificing valuable laboratory space.**

## Benefits of the Access 2 immunoassay system

### **Maximized productivity in a minimal footprint**

- › Throughput of up to 100 tests per hour
- › One of the smallest footprints in its class
- › Loading capacity of 60 samples at one time
- › Random access for running any combination of samples at any time
- › Rapid kinetics for improved turnaround time
- › STAT and automatic reflex testing

### **Reliable uptime with fewer emergency service calls**

- › Sample probe obstruction detection
- › Remote monitoring and instrument diagnostics capability through PROService

### **Established performance for confidence in patient results**

- › Utilizes Proven chemiluminescent and magnetic particle separation technology
- › Unique, assay-specific protocols for optimal assay performance

### **Lower operating costs**

- › Onboard refrigerated storage of up to 24 assays
- › 50 tests per pack\*
- › No dedicated water system required

### **Intuitive interface and improved ergonomics simplify operation**

- › Customer-friendly user interface with color-coded, icon-driven touch screen
- › Internal barcode reader for automatic identification of reagent packs upon loading, minimizing operator loading errors
- › Ergonomic workstation with flexible movement for comfort for system operators
- › Status bar indicator light for easy visibility of system status from across the laboratory

## Standardization simplifies and improves operational efficiencies

All Access immunoassay systems feature standardized test menus, assay protocols, instrument processes and reference ranges. Together, these features minimize training requirements and streamline inventory management tasks for labor-cost savings, while providing consistent results across the healthcare delivery network.

## Comprehensive training and top-tier support increases productivity

Beckman Coulter is committed to ensuring its diagnostic testing solutions meet the highest standards and integrate seamlessly into different laboratory environments. The Access 2 comes complete with the training and support needed to get it up and running as quickly as possible.

- › Training provided (advanced operator training also available)
- › Online maintenance help
- › Dedicated customer service and support to meet current and future needs



## The Beckman Coulter Diagnostics Difference

Beckman Coulter has invested in a synergistic portfolio of products and services for modern healthcare enterprises that are scalable to fit laboratories of any size. This portfolio includes more than clinical diagnostic products and services; it also includes clinical IT solutions, the Danaher Business System, and our top-tier service and uptime solutions. These four elements work together in synergy to help our customers achieve laboratory excellence. We call it “the Beckman Coulter Diagnostics Difference.”

DANAHER BUSINESS SYSTEM

CLINICAL INFORMATION  
MANAGEMENT TOOLS



COMPREHENSIVE, HIGH-QUALITY  
CLINICAL DIAGNOSTIC SOLUTIONS

TOP-TIER SERVICE AND UPTIME

# ACCESS IMMUNOASSAY MENU

## Adrenal/Pituitary

- > Cortisol (Serum, Plasma, Urine)

## Anemia

- > Active B12\*
- > EPO
- > Ferritin
- > Folate
- > Folate (RBC)
- > Intrinsic Factor Ab
- > sTfR
- > Vitamin B12

## Bone metabolism

- > 25(OH) Vitamin D Total
- > Intact PTH  
(Routine or Intraoperative)
- > Ostase
- > Ultrasensitive hGH

## Cardiac

- > AccuTnI+3 (Troponin I)
- > CK-MB
- > Digoxin
- > hsTnI\*
- > Myoglobin
- > Triage BNP\*\*

## Diabetes

- > Ultrasensitive Insulin

## Infectious disease

- > Rubella IgG
- > Toxo IgG
- > Toxo IgM II

## Reproductive

- > AFP (ONTD)
- > AMH
- > DHEA-S
- > hFSH
- > hLH
- > Inhibin A
- > Progesterone
- > Prolactin
- > Sensitive Estradiol
- > SHBG
- > Testosterone
- > Total  $\beta$ hCG (5<sup>th</sup> IS)
- > Unconjugated Estriol

## Thyroid

- > T3 (Free and Total)
- > T4 (Free and Total)
- > Thyroglobulin
- > Thyroglobulin Ab II
- > Thyroid Uptake
- > TPO Ab
- > TSH (3<sup>rd</sup> IS)

## Tumor markers

- > AFP
- > CEA
- > BR Monitor  
(CA 15-3 Antigen)
- > GI Monitor  
(CA 19-9 Antigen)
- > OV Monitor  
(CA 125 Antigen)
- > Prostate Health  
Hybritech PSA  
Hybritech free PSA  
Hybritech p2PSA  
Prostate Health Index  
(*phi*)\*\*\*

## Research use only

- > PAPP-A (RUO)
- > IL-6 (RUO)

## Access immunoassay menu


Beckman Coulter offers an extensive immunoassay diagnostic testing menu to fulfill patient-care needs.

- > More than 50 diagnostic immunoassays available
- > Multiple assays available for different diagnostic specialties

## Easy-to-use, easy-to-store reagent kits

The Access 2 uses liquid, ready-to-use reagent packs that contain everything needed to run an assay in one convenient container.

- > Slim-line design, facilitating easy pack loading
- > Self-sealing elastomer closure
- > No caps, preventing potential probe crashes
- > Easy-to-store kit boxes maximize refrigerator space
- > Typical shelf life of 12 months

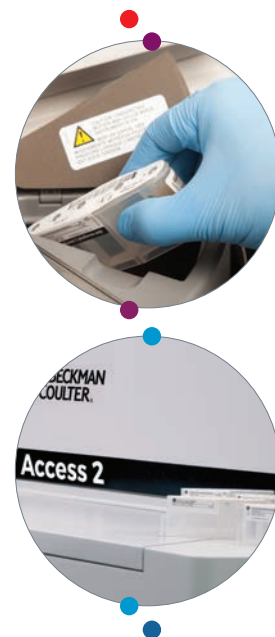

+ TSH (3<sup>rd</sup> IS) is available in 100-test pack.

\* In development.

\*\* Available exclusively from Quidel or its authorized distributors for use on Beckman Coulter immunoassay systems.

\*\*\* *phi* is a combination of Hybritech PSA, free PSA and p2PSA.

Brochure is intended for the United States.

Contact your local Beckman Coulter sales representative for availability.

© 2018 Beckman Coulter, Inc. All rights reserved. Beckman Coulter, the stylized logo and the Beckman Coulter product and service marks mentioned herein are trademarks or registered trademarks of Beckman Coulter, Inc. in the United States and other countries.

For Beckman Coulter's worldwide office locations and phone numbers, please visit [www.beckmancoulter.com/contact](http://www.beckmancoulter.com/contact)